

Portable suction pump

Features


- Oilless and maintenance free vacuum pump
- 2-liter waste collector with overflow protection
- On / off foot switch for vacuum pump
- Available in 110 or 220 Volts

Specifications

- Weight: 33 lbs (15 kg) for 2.6 CFM
41 lbs (18.5 kg) for 4.6 CFM
- Noise level: 44 Db (for both motor)
- Dimensions: 19.70" x 12.00" x 18.00" (50 x 30.5 x 45.7 cm)
- Air displacement vacuum: 2.6 CFM (73.6 L/min) or 4.6 CFM (127.4 L/min)

Options

- S.E. valve with asepsis tubing | motor 2.6 CFM (SAF-606-2.6)
- HVE valve with asepsis tubing | motor 2.6 CFM (SAF-606-2.6)
- S.E. **and** HVE valve with asepsis tubing | motor 4.6 CFM (SAF-606-4.6)

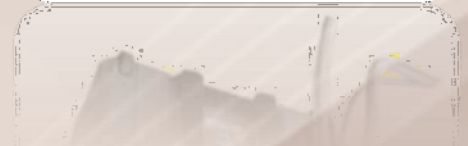

Optional handpiece control for portable suction

Features

- Air / water syringe with quick change tip
- Midwest 4-hole connection for use of a low speed handpiece (dry)
- Midwest 4-hole connection for use of a high speed handpiece (wet)
- Control of the handpieces through a rheostat foot pedal
- 1 litre bottle for the clean water system
- All tubing are easy to clean, and instruments are autoclavable
- 1/4" Male air quick disconnect with 7' tubing for add-on compressor
- External quick disconnect water for auxiliaries such as Cavitron
- Handpieces water flow control adjustment

Specifications

- Weight: 40 lbs (18.15 kg) for 2.6 CFM | 48 lbs (21.75 kg) for 4.6 CFM
- Dimensions: Fit inside de portable suction pump



One-year warranty on manufacturing defects